

SOMNI  
SCIENTIFIC



## VETERINARY EXPOSURE PREVENTION SYSTEM

Somni Scientific Veterinary Exposure Prevention System is a complete wall unit evacuation system. Powered by a vacuum fan unit it will effectively scavenge waste anesthetic gas from a single or multiple anesthesia systems. Utilizing 1.5" PVC it can be connected to wall plates or ceiling plates. It works in conjunction with the Somni Interface Valve to provide the ultimate control/draw from your absorber and/or non re-breather.

### ► SPECIFICATIONS

<b>WEIGHT: 7.5 lbs.</b>	<b>POWER: 110v</b>
<b>HEIGHT: 9.74 in.</b>	<b>Db: 20</b>
<b>DEPTH: 6.25in.</b>	<b>INLET: 1.5"</b>
<b>WIDTH: 7.74 in.</b>	<b>OUTLET: 2"</b>




### WALL UNIT FEATURES

- Lighted On/Off switch
- Hour Meter
- Easy external mounting points
- Ultra quiet operation



### INTERFACE FEATURES

- Manufactured from high quality, durable materials.
- Standard 19mm inlet and outlets
- Multiple mounting options
- Dual inlets for re-breathing and non re-breathing
- Adjustable air intake for balancing

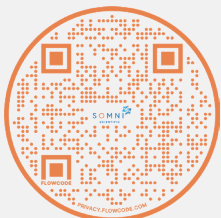

1-877-637-3625

info@somniscientific.com

www.somniscientific.com

Proudly machined, assembled, and calibrated in the USA