

SOMNI
SCIENTIFIC



AMD 3

ACCURATE MULTIPLE DELIVERY

Somni Scientific's AMD (Accurate Multiple Delivery)

The Somni AMD-3 accurate multiple delivery is a flowmeter controlled diverter system that is adaptable to any standard Inhalant Anesthesia System. It can be added to Somni's TT-2 or any other non-pressurized system.

The AMD-3 gives the user the ability to deliver accurate flowrates to 3 stations simultaneously. No more guessing if the flow is being delivered equally to multiple devices at once. With its included pressure relief valve the user does not have to worry about pressure building up in the vaporizer. 3 fresh gas common outlets supplied with color coded tubing makes this system very user friendly.

► SPECIFICATIONS

Weight: 11 lbs / 5 kg

Height: 7 in / 177.8 mm

Depth: 8 in / 203.2 mm

Width: 12 in / 304.8 mm




FEATURES:

- 3 Flowmeter Controlled Stations with Toggle Switches
- Color Coded Tubing to each station
- Common Gas outlets for each station
- Small footprint
- Adjustable Flowrates
- Stacking design saves valuable workbench space



INCLUDES:

- Pressure Relief Valve
- 3 each, 6' color coded tubing with adapters

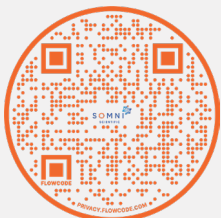

1-877-637-3625

info@somniscientific.com

www.somniscientific.com

Proudly machined, assembled, and calibrated in the USA