

SOMNI
SCIENTIFIC





NRB SYSTEM

NRB System (Non Re-breathing System)

Somni Scientific's Passive style Non Re-Breathing Mouse or Rat Nose Cone System was developed to help minimize exposure to users while also adding key features to simplify the process for the technician.

A magnetic base has been added to the nose cone eliminating the need to tape down the device. Magnetic animal positioners are also included to help the user keep the animal in the ideal position during their procedure. 10mm thick latex material provides a proper seal around the rodent's nose, keeping the gases inside the nose cone and minimizing waste gas exposure in the workplace.

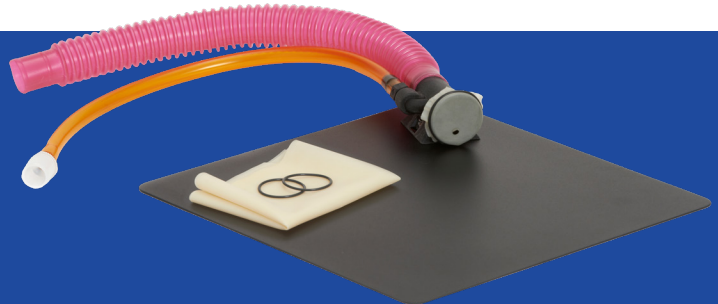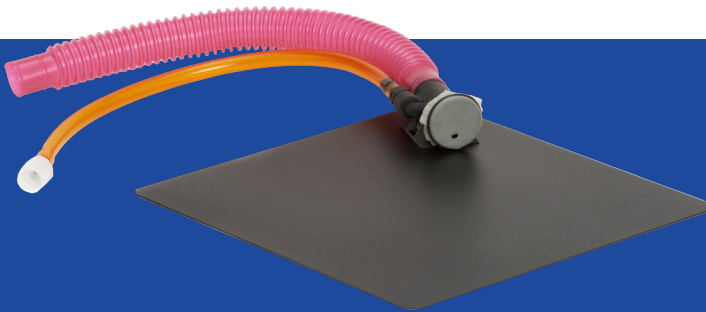

FEATURES:

- Magnetic Base
- Magnetic Animal Positioners
- Low Flow Design
- 19mm or 15mm for Waste Gas Scavenging
- Color Coded Tubing



INCLUDES:

- 12x12 Coated Steel Work Base
- 2'x 6" 10mm Latex material
- 5 O-Rings to hold Latex
- 2 Animal Positioners



1-877-637-3625

info@somniscientific.com

www.somniscientific.com

Proudly machined, assembled, and calibrated in the USA